



POSTERIOR CROWN

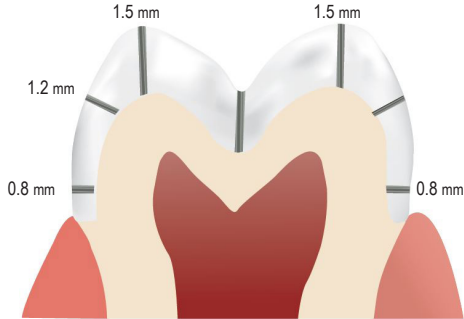

BUTT MARGIN

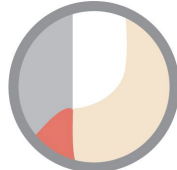

CHAMFER




CEMENTATION

Conventionally cement following the manufacturer's instructions for use with Zirconia restorations.


CONNECTOR CROSS SECTIONS

Anterior	$\geq 12\text{mm}^2$
Posterior	$\geq 12\text{mm}^2$

ANTERIOR CROWN

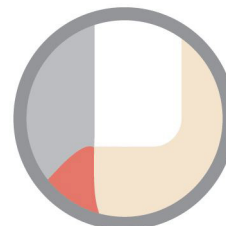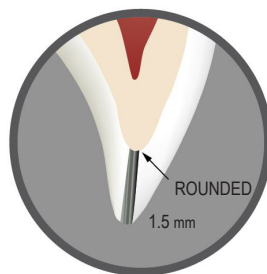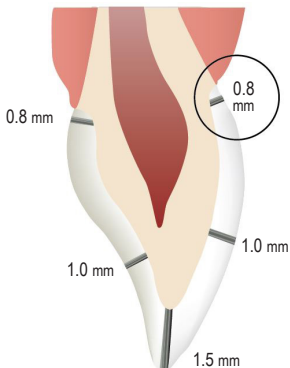

BUTT MARGIN

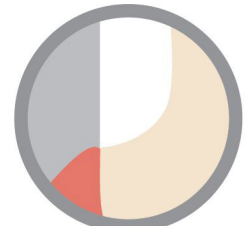

CHAMFER

INDICATIONS FOR USE

		Zirconia
Anterior	Crown	✓
	Bridge (up to 3 units)	✓

		Zirconia
Posterior	Crown	✓
	Bridge (up to 3 units)	✓

CHAIRSIDE ADJUSTMENT

If chairside adjustments are needed:

- Fine diamond tools are recommended using a water-cooled, high speed hand piece.
- Use light pressure to avoid overheating and micro-cracking of the sintered body.
- Polish with diamond impregnated rubber polishers; maintain minimum wall thickness and connector cross sections.